

IMS Series Automatic Flake Ice Maker

Instruction Manual

Thank you for selecting our IMS Series Automatic Flake Ice-Maker!

To ensure the proper operation of this ice maker, please read the instruction

carefully before use and please keep it for inspection in case troubles occur.

CONTENT

FEATURES

NOTES BEFORE OPERATION

STRUCTURE

OPERATION PROCEDURES

MAINTENANCE

TECHNICAL PARAMETERS

TROUBLE AND SHOOTING

CIRCUIT DIAGRAM

FEATURES:

Our IMS-series automatic flake ice makers are designed in cubic shape with the curved front panel and curved door (or stainless steel flat front panel and flat door) for nice appearance. Many advantages of the ice maker are available, such as ice making quickly, great ice making capacity, nice ice shape and ice dropping quickly. Tap water inlet → water incoming → ice making → crush ice → ice dropping → ice storing, all these processes are automatically controlled for ice making continuously. In case of water shortage or ice full in cabinet, the indicator on the operation board will light on accordingly and the ice maker will stop operation automatically. The ice storage cabinet is PU foamed, so it is insulated well and it could prevent the ice made from melting. Also, the remaining water after ice made each time is used recycling. So not only the water is saved, but also the ice making efficiency is enhanced and the energy consumption is reduced, and the customer's ice making cost is also reduced much. Either purified water or the tap water could be used for ice making.

NOTES BEFORE OPERATION:

1. The incline angle of the cabinet could not be over 45° during transportation. Don't make ice maker upside-down in case the compressor or refrigerating system troubles occur.

2. The flake ice maker should be placed on horizontal and structurally sound ground and away from heat source and corrosive gas. At least about 150mm room should be left around the ice maker to ensure the fine ventilation.

3. Before using the flake ice maker for the first time, you should wait 12 hours after positioning it in a proper place.

4. Voltage range:187~242V(220V)

100~130V(110V)

5. The rated capacity of wire should be over 6A,the section area of the wire is 0.75mm².The wire could be consisted by single ply or multiplies.

The rated current of fuse is 5A.


6. The separate three-pole outlet should be used and the grounding must be good.

7.If the supply cord is damaged,it must be replaced by the manufacture,its dervice agent or similarly qualified

persons in order to avoid a hazard.

8. This appliance is not intended for use by children or other persons without assistance or supervision if their physical, sensory or mental capabilities prevent them from using it safely. Children should be supervised to ensure that they do not play with the appliance.

STRUCTURE:

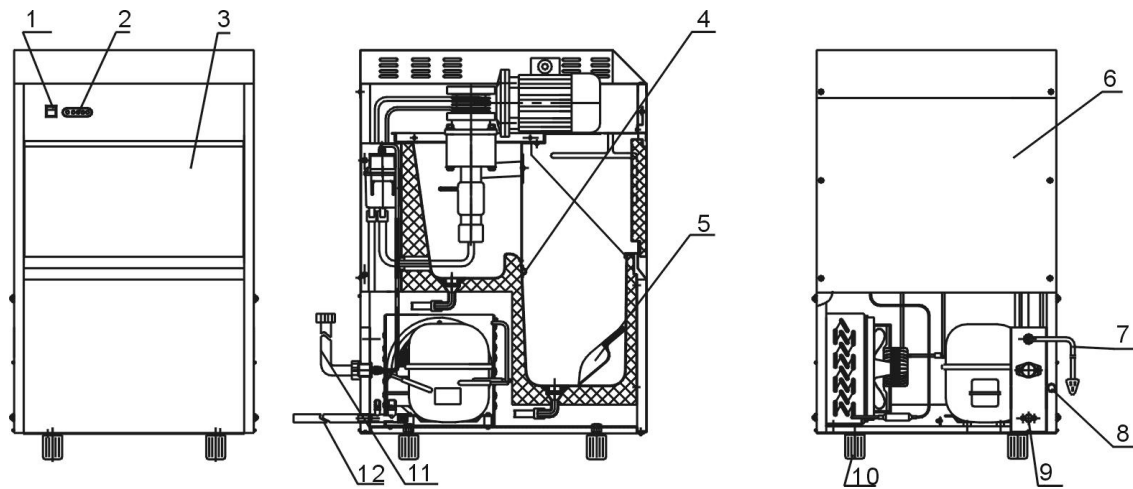

1. Switch
2. Operation Board
3. Door
4. Temperature Sensor
5. Ice Shovel
6. Back Panel
7. Power Cable
8. Screw For Drainage

9. Water Outlet Tube
10. Foot Screw
11. Water Inlet Tube
12. Drainage Tube

OPERATION PROCEDURES:

1. Unpack the package and take out the attached documents and spare parts such as the water inlet and outlet tubes, ice shovel and gaskets, etc.
2. Position the ice maker on a well ventilated place and leave at least 150mm room between the ice maker and the wall. Be sure the ice maker is positioned levelly and away from the heat source.
3. Make one end of the $\phi 12$ flexible plastic corrugated pipe (supplied with the maker) connected with the water outlet tube on back of the maker, another pipe end should be placed in a container for surplus water (prepared by user himself) or in the sewer.
4. Make one end of the water inlet tube (supplied with the maker) connected with the 3/4" screw type connector of the water tap so as to get the drinkable tap water. The

water pressure of water supply pipe is 1.5 ~ 3 bar. The other end of the water inlet tube should be connected with the screw type connector of water valve on the back of the ice maker. Be sure to put gaskets (supplied with the maker) in both ends of the water inlet tube before connecting. If it is an ice maker for purified water, a $\phi 12$ corrugated pipe (supplied with the maker) should be used. Connect one end of the pipe with the purified water supplier (desk-top water dispenser) and another end with water valve connector on the back of the ice maker. The desk-top water dispenser should be placed on the ice maker at the back position.

5. Make power connection and press on the start button on the operation board for about 4 seconds till the run indicator and the big ice indicator (or small ice indicator) light on. Then the maker begins to work. All the procedures, water inlet → ice making → crush ice → ice dropping → ice storing, are controlled automatically for ice making continuously. If ice is full in the ice storage container, the ice full indicator on the operation board will light on and the ice maker will stop automatically. If there is a water shortage or something wrong with the

water supply system, the Water shortage indicator will light on and the ice maker will stop operation automatically.

MAINTENANCE:

1. If the compressor stops for any reasons like water shortage, too much ice, power off, etc., don't restart it right away.

You can restart it about 5 minutes later.

2. Check regularly the connectors of the water inlet and outlet tubes and drain the little surplus water that may leak.

3. If the ice make will be unused for a long period, please screw off the plastic drainage connector screw on the back of the ice maker so as to drain the surplus water in the water groove. Then screw it on. Wipe the inner liner of the ice storage container with a clean rag.

4. When plugging or unplugging, the plug should be held by hand and the wires should not be dragged heavily.

TENCHNICAL PARAMETERS

| Model | Ice making capacity (Kg/24h) | Power input (W) | Refrigerant R134a (g) | Dimensions (mm) |
|----------|------------------------------|-----------------|-----------------------|-----------------|
| IMS-20 | ≥ 20 | 180 | 100 | 300 × 493 × 547 |
| IMS-40 | ≥ 40 | 280 | 130 | 380 × 543 × 722 |
| IMS-50 | ≥ 50 | 280 | 130 | 380 × 543 × 722 |
| IMS-70 | ≥ 70 | 420 | 140 | 548 × 611 × 883 |
| IMS-85 | ≥ 85 | 420 | 140 | 548 × 611 × 883 |
| IMS-100 | ≥ 100 | 420 | 150 | 548 × 611 × 883 |
| IMS-150S | ≥ 150 | 685 | 250 | 548 × 611 × 883 |
| IMS-150B | ≥ 150 | 685 | 250 | 635 × 611 × 945 |

The specifications are subject to be changed without notice.
Please check the nameplate for sure.

***TROUBLE AND SHOOTING
(FOR REFERENCE OF USERS AND TECHNICIANS)***

| TROUBLE | CAUSE | REMEDY |
|--|---|--|
| The flake ice maker does not work. | <ol style="list-style-type: none"> 1. The voltage is lower than the lower limitation 2. The ambient temperature is too low (Lower than 10°C) | <ol style="list-style-type: none"> 1. Stop the ice maker and restart it until the voltage is normal. 2. Try again when the ambient temperature higher than 10°C. |
| Water shortage indicator light on. | <ol style="list-style-type: none"> 1. Water shortage from the water supply system. 2. Water valve damage. 3. Tap water pressure too low. | <ol style="list-style-type: none"> 1. Check the water supply. If it is OK, restart the maker. 2. Check the water valve. 3. Be sure the tap water pressure is higher than 1.5 bar. |
| Compressor does not work. | <ol style="list-style-type: none"> 1. Water shortage. 2. Ice full. | <ol style="list-style-type: none"> 1. Check the water supply system (water valve and water supply tube). 2. Take off some ice. |
| Compressor works, but no ice is made. | <ol style="list-style-type: none"> 1. Refrigerant leakage. 2. The refrigerating system is blocked. | <ol style="list-style-type: none"> 1. Recharge the refrigerant and check the leakage again. 2. Charging N₂ to check whether the system is blocked and replace the filter. 3. Check whether the fan is running. |
| Water shortage indicator and ice full indicator are glittering at the same time. | If water inlet jam or running failure because of low environmental temperature. | Stop the machine until the temperature reach 10°C or repair it by authorized service people. |

ICIRCUIT DIAGRAM:

